IR Thermometer

Instruction Manual

**DIFRANCE**<sup>®</sup>

CTG100

# English

## Warning and Precautions:

Keep out of reach of children under 12 years.

Never use the thermometer for purposes other than those for which it has been intended. Please follow the general safety precautions when using on children.

Never immerse the thermometer into water or other liquids (not waterproof). For cleaning and disinfecting please follow the instructions in the <<Care and cleaning>> section.

Do not store this thermometer in temperature extremes below -25 °C or over 55 °C (below -13 °F or over 131 °F) or in excessive humidity (above 95% non-condensing relative humidity).


If thermometer is stored in a location that is cooler or warmer than where it will be used, allow it to acclimate to the room temperature for 10 minutes before taking a measurement.

Do not use the thermometer if there are signs of damage on the scanner or on the thermometer itself. If damaged, do not attempt to repair the product.

Never insert a sharp object into the scanner area or any other open surface on the thermometer.

This thermometer consists of high-quality precision parts. Do not drop the instrument. Protect it from severe impact and shock. Do not twist the instrument or the measuring sensor.

This thermometer is intended for household use only.
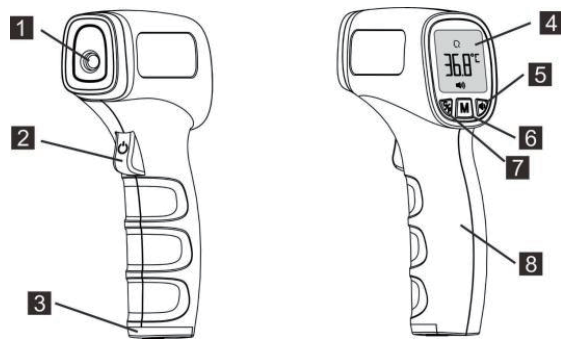
## Location of controls:

1. IR Sensor
2. Scan button/ Power on button
3. Battery Compartment
4. LCD display
5. Unit sound on/ Mute button
6. Mode button / Switch of Body temperature or surface temperature
7. Celsius °C/ Fahrenheit °F switch button
8. Handgrip

## Measurement in under 2 seconds

The innovative infrared technology allows measurement of temperature in under 2 seconds whether you use the body or surface temperature measurement option.

## Accurate and reliable

Due to the unique technology, the IR forehead thermometer captures the heat naturally given off by the forehead to calculate body temperature value. Get the same professional accuracy whether touching your forehead or holding up to 2.5 cm (1 inch) away.

## Easy to use

The No touch + forehead thermometer is non-invasive. A measurement can be taken even while a child is sleeping.

The No touch + forehead thermometer is safer to use on a child than a rectal thermometer and is easier to use than other methods.

## Safe and hygienic

No touch option helps minimize spreading of germs. Completely safe for use on children and adults.

## Power on

Open the battery compartment and insert 2pcs x AA size (1.5V LR6) batteries (not included) into the battery compartment.

Press Power on button to start.

## Scan Temperature

Simply press "SCAN" button to scan the temperature

## Body temperature and Surface temperature Switch

Simply press MODE button to switch

## More functions:

Memory saved

The unit will save the temperature automatically every time. There are totally 32 memories can be saved, the unit will save the last 32 memories.

## F-1

Power on the unit. Press and hold SET button for 3 seconds to enter parameter settings, when F1 is appeared, press °C/ °F and Sound on/Mute button to review the saved memories.

## F-2

While in F-1 mode, press MODE button again to go to F2. Press the °C/°F or the Sound on/mute button. The fever alert threshold is displayed. The threshold value will be changed by 0.1°C/°F every time when the °C/°F button or sound on/mute button is pressed. The tunable range is 35.0°C–42.0°C (95.0°F–107.6°F). The default fever alert threshold is 37.6°C. If the body temperature is measured higher than 37.5°C, the LCD display will turn RED light. The RED light will be appeared if the measured body temperature is higher than the set fever alert threshold.

## Screen display showing:

HI

In Object mode, a temperature reading of more than 100°C

In Forehead mode, a temperature reading of more than 42.2°C (108.0°F)

Lo

In Object mode, a temperature reading of less than 0°C (32.0°F)

In Forehead mode, a temperature reading of less than 34.9°C (94.8°F)

## To RESET the unit

1. Power on the unit. Press and hold the °C/°F button and Sound on/mute button at the same time for 3 seconds.
2. "dEL" will show on LCD display. Then press the Scan button.
3. "End" will then show on LCD display. All data were cleared.

## Specifications

Product Name	Infrared Thermometer
Product Model	CTG 100
Power Supply Mode	Battery operated
Operating Voltage	DC 3V
Battery Model	2pcs x AA sized batteries
Operating Mode	Continuous operating
Display	Segment LCD
Measure time	About 1 second
Latency Time	About 1 second
Emissivity	0.95
Measuring Distance	1 to 3 cm
Measuring Range	Forehead: 34.9°C–42.2°C (94.8°F–108.0°F) Object: 0.0°C–100.0°C (32.0°F–212.0°F)
Accuracy (Laboratory)	±0.4°F/±0.2°C from 94.8°F to 108.0°F (34.9°C to 42.2°C) ±0.6°F/±0.3°C, Outside the range of 94.8°F to 108.0°F (34.9°C to 42.2°C)
LCD backlight	Green: Body temperature is below 37.6°C (default), or below the set fever alert threshold RED: Body temperature is higher than 37.5°C (default), or higher than the set fever alert threshold White: In Surface temperature mode
Memory	32 memories
Automatic Power-off	The thermometer automatically powers off if it is not used in 10 seconds.
Dimensions (mm)	150 x 92 x 46
Operating Environment	Temperature: 10°C–40°C (50°F–104°F) Humidity: 15%–95% RH, non-condensing
Storage and Transportation	Temperature: -20°C to 55°C (-4°F–131°F) Humidity: 0- 95% RH, non-condensing

Use of this thermometer is not intended as a substitute for consultation with your physician.

Temperature elevation may signal a serious illness, especially in neonates and infants, or in adults who are old, frail, or have a weakened immune system. Please seek professional advice immediately when a temperature elevation occurs on persons that are:

- Neonates and infants under 3 months (Consult your physician immediately if the temperature exceeds 37.4 °C [99.4 °F])
- Individuals over 60 years of age (Fever may be blunted or absent in older individuals)
- Individuals having diabetes mellitus or a weakened immune system (e.g., HIV positive, cancer chemotherapy, chronic immunosuppressant treatment, splenectomy)
- Individuals who are bedridden (e.g., nursing home patient, stroke, chronic illness, paraplegia, quadriplegia, surgical recovery)
- A transplant recipient (e.g., liver, heart, lung, kidney).

This thermometer is not intended for pre-term babies or small-for-gestational- age babies. Pre-term is defined as babies with a calendar age of less than 37 weeks. Small-for-gestational-age is defined as a baby born at 37 weeks or later, with a weight below the 10th percentile for babies of the same gestational age.

This thermometer is not intended to interpret hypothermic temperatures. If the device displays a temperature of 36.4 °C (97.5 °F) or less and the individual is exhibiting atypical symptoms or behaviors, contact your physician or health care professional.

DO NOT allow children to take their temperatures unattended.

Please consult your physician if you see symptoms such as unexplained irritability, vomiting, diarrhea, dehydration, changes in appetite or activity, seizure, muscle pain, shivering, stiff neck, pain when urinating, etc., even in the absence of fever.

Even in the absence of fever, those who exhibit a normal temperature may still need to receive medical attention. People who are on antibiotics, analgesics, or antipyretics should not be assessed solely on temperature readings to determine the severity of their illness.

DO NOT modify this equipment without the authorization of the manufacturer.

Please notice - All products are subject to change without any notice. We take reservations for errors and omissions in the manual.

ALL RIGHTS RESERVED, COPYRIGHT INTER SALES A/S

# DIFRANCE<sup>®</sup>

CE

Electric and electronic equipment contains materials, components and substances that can be hazardous to your health and the environment, if the waste material (discarded electric and electronic equipment) is not handled correctly.

Electric and electronic equipment is marked with the crossed out trash can symbol, seen above. This symbol signifies that electric and electronic equipment should not be disposed of with other household waste, but should be disposed of separately.

All cities have established collection points, where electric and electronic equipment can either be submitted free of charge at recycling stations and other collection sites, or be collected from the households. Additional information is available at the technical department of your city.

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ALLE RECHTE VORBEHALTEN, URHEBERRECHTE INTER SALES A/S
Inter Sales A/S